
S950 / S900 Air Scaler Tips




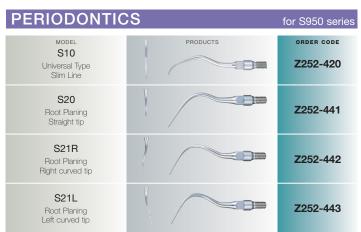

*Note: S10 / S20 / S21R / S21L must be used with the NSK S950 series.

ENDODONTIC	S V-Tip	for S950 series
MODEL S75 V-Tip holder	PRODUCTS	ORDER CODE Z252-035
V-U15 V-file irrigation 23mm #15 ISO15	(Pack of 6)	Y900-177
V-U20 V-file irrigation 23mm #20 ISO20	(Pack of 6)	Y900-178
V-U25 V-file irrigation 21mm #25 ISO25	(Pack of 6)	Y900-179
V-U30 V-file irrigation 22mm #30 ISO30	(Pack of 6)	Y900-180
V-U35 V-file irrigation 23mm #35 ISO35	(Pack of 6)	Y900-181
V-N50 Ni-Ti irrigation 17mm	(Pack of 3)	Y900-182
V-N51 Ni-Ti irrigation 21mm	(Pack of 3)	Y900-183

*Note : Endodontic V-Tips must be used with the NSK S950 series.

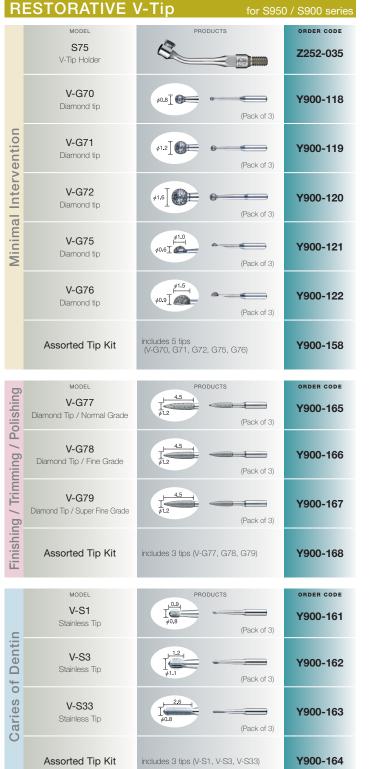

KaVo and MULTiflexLUX* are registered trademarks of Kaltenbach & Voigt GmbH & Co., Germany. Sirona* is a registered trademark of Sirona Dental Systems GmbH, Germany. W&H* and Roto Quick* are registered trademarks of W&H Dentalwerk Bürmoos GmbH, Austria. Bien-Air* and Unifix* are registered trademarks of Bien-Air Dental S.A., Switzerland.





All NSK handpieces are built to have excellent durability against repeated autoclaving at temperatures of 135°C.



NAKANISHI INC. 700 Shimohinata, Kanuma, Tochigi 322-8666, Japan TEL:+81 (0)289-64-3380 FAX:+81 (0)289-62-5636 www.nsk-inc.com

TOKYO OFFICE 6F MY Building, Ueno 1-13-3, Taito-ku, Tokyo 110-0005, Japan TEL:+81 (0)3-3835-2892 FAX:+81 (0)3-3835-2856

NSK Europe Westerhachstraße 58 D-604

Westerbachstraße 58 D-60489 Frankfurt, Germany TEL:+49 (0)69 74 22 99 0 FAX:+49 (0)69 74 22 99 29 www.nsk-europe.de



Dubai Airport Free Zone, Room 2W-310, 3F West Wing, Dubai, UAE TEL:+971-4-299-6644 FAX:+971-4-299-6664



TI-MAX S950 / S900 SERIES FROM NSK

MULTI-PURPOSE AIR SCALER



compact, dynamic and now with the unique power ring pre-set power adjustment





	APPLICATIONS	TIPS	NSK POWER RING SETTING LEVEL 1 LEVEL 2 LEVEL 3
Sc	aling	S1, S2, S3	
Pe	riodontics	S10, S20, S21R, S21L	
En	dodontics	S75 + V tips (U-file, Ni-Ti)	
Restorative	Minimal Intervention	S75 + V tips (V-G70, 71, 72, 75, 76)	
	Finishing, Triming, Polishing	S75 + V tips (V-G77, 78, 79)	
	For Caries of Dentin	S75 + V tips (V-S1, S3, S33)	



Specific Power Adjustment

Power Ring allows adjustment to 3 specific power settings - from the gentle minimum setting for delicate periodontics procedures to maximum power for general scaling.



Optic Illumination

S950 emits a bright circular light providing clear operator vision, without shadow, direct to the operating field.



Solid Titanium Body

Constructed from Solid Titanium, the S950 body is light weight, very well balanced and extremely durable. The actual characteristic of titanium material provides a firm, tactile grip, and offers a high degree of operator control and comfort.



Tip Wrench with Torque Limiter

The tip wrench has an integrated torque limiter. The wrench completely eliminates tip breakages by preventing over tightening of tips during tip exchange.



The Extensive Tip Range

The wide range of tips provides applications for a variety of clinical procedures including Periodontics, Endodontics, Scaling and MI. (Refer page 4)



SCALING

ENDODONTICS

S950 / S900 Air Scalers

Ti-Max S950 Series with Unique Power Ring

Ti-Max S950 is multi-purpose air scaler which, in combination with the extensive tip range, provides a wide variety of clinical applications. 3 power levels can be selected via the unique NSK Power Ring. A single S950 air scaler has the capability for Periodontics, Endodontics, Scaling and MI procedures. The Circular Optic illumination increases intraoral vision to the operating field. S950 has a Solid Titanium body which is durable and lightweight. The scaler offers excellent adjustable power, but low noise. Ti-Max S950 series directly connects to the same NSK Phatelus coupling as for high-speed handpieces. Available in optic and non-optic models.



Ti-Max S950KL / S950SL / S950BL / S950WL with Unique Power Ring

Ti-Max S950KL / S950SL / S950BL / S950WL series is titanium body air scaler with circular light. There power levels can be selected via the unique NSK Power Ring. Ti-Max S950KL / S950SL / S950BL / S950WL directly and easily connect to other brands coupling as for high speed handpieces.



Ti-Max S900 Series without Power Ring

Ti-Max S900 is follows the traditional NSK Air Scaler features, without the Power Ring, but upgraded to a Solid Titanium body. S900 offers excellent power, low noise and connects directly to the NSK Phatelus coupling. The circular light emission clearly illuminates the operational area facilitating the possibility of faster treatment with reduced operator stress. Each scaler is delivered complete with a starter set of three tips including universal, sickle, and Perio. Available in optic and non-optic models.



Solid Titanium body / Circular optic (S900L) / Frequency of oscillation: 6,100~6,600 Hz



• Weight: 63a • includes 3 tips

*Both Ti-Max S950 and S900 series includes 3 tips (S1, S2, S3), CR-2 Tip Wrench and Tip Cover L *Only external body component is Titanium.